

Copure® Multi-Functional Purification Column

# Fast Mycotoxin Determination Solution

Copure® Multi-Functional Purification Column contains multiple adsorption matrices to adsorb quickly and selectively of lipids, organic acids, proteins, pigments and other impurities. The matrices do not absorb the target to realize fast purification.

## Features


- Complete purification simply and fast in 30s
- Stored at room temperature for more than 18 months
- Recovery rate  $\geq 90\%$ , RSD  $\leq 5\%$



# Copure® Multi-Functional Purification Column Complete purification within 30s

## Procedures

Add sample to the tube, then push the purification column to the bottom of the tube. The purified sample passes through the filler to the top. After dry and re-dissolve, the target could be tested.

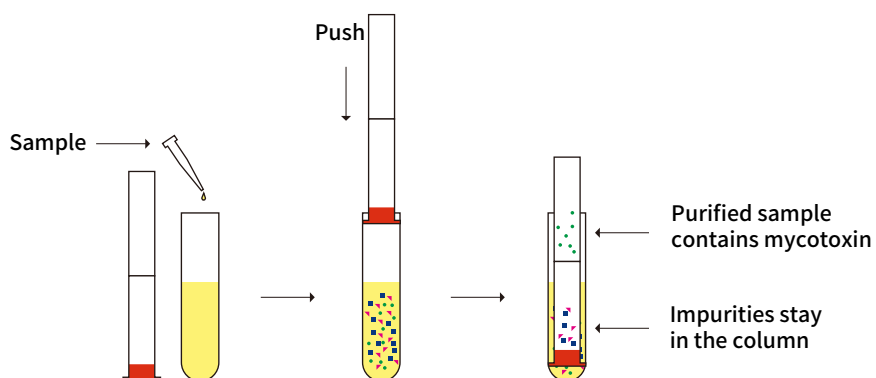

## Ordering Information:

Cat. #	Description	Application	Qty.
COAF228	Copure® 228 Multi-Functional Purification Column	Patulin, Aflatoxin B1 B2 G1 G2	25 pcs/box
COAF226	Copure® 226 Multi-Functional Purification Column	Zearalenone, Aflatoxin B1 B2 G1 G2	25 pcs/box
COAF224	Copure® 224 Multi-Functional Purification Column	Zearalenone	25 pcs/box
COAF223	Copure® 223 Multi-Functional Purification Column	Aflatoxin M1 M2	25 pcs/box
COAF230	Copure® 230 Multi-Functional Purification Column	Deoxynivalenol	25 pcs/box
COAF229	Copure® 229 Multi-Functional Purification Column	Ochratoxin A	25 pcs/box
COAF302	Copure® 302 Multi-Functional Purification Column	Multiple functions	25 pcs/box

**BIOCOMMA LIMITED**

**Add.:** Ground FL, Bldg. 12, Zhonghaixin Innovation Industrial Center, 12 Ganli 6th Rd., Jihua St., Longgang Dist., Shenzhen, Guangdong, 518114 P.R. China

**Tel.:** +86 (0755)-25431879

**Web:** www.biocomma.com

**Email:** info@biocomma.com